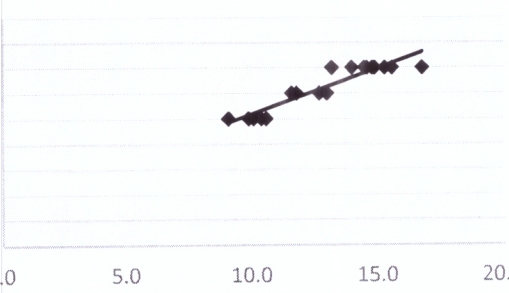


Batch Evaluation Report

Name of the Test Evaluated	:	Colour Scale for Hemoglobin (Hemocheck)
Name of Manufacturer	:	Allied Health Sciences Pvt Ltd F-95, Okhla Industrial Estate Phase III New Delhi - 110020
Batch Number	:	27//9999
Manufacturing Month/Year	:	11/2023
Expiry	:	12/2029
Evaluated During	:	November 2023
Evaluation Protocol	:	➤ Comparison of twenty specimens with Mindray 5 Part BC 6200 Hemoglobin values
Acceptability Limit	:	$r > 0.95$, $p > 0.05$

Batch Evaluation Report

Results (contd)	Correlation Curve $y = 0.7195x + 3.2506$ $r^2 = 0.8745$ $r = 0.94$ 
Conclusion	: Correlation of Hemocheck HB with Automated Cell Counter Mindray BC 6200 was found to be good ($r = 0.94$, $p > 0.05$)

Prepared by:



Mr. Rahul Chavan
Junior Research Officer

Reviewed by:



Dr. Rashmi D Khadapkar
Sr. DGM – R&D